


Alpha - 660 Dual-Disc Automatic Grinder and Polisher

GRINDING/POLISHING

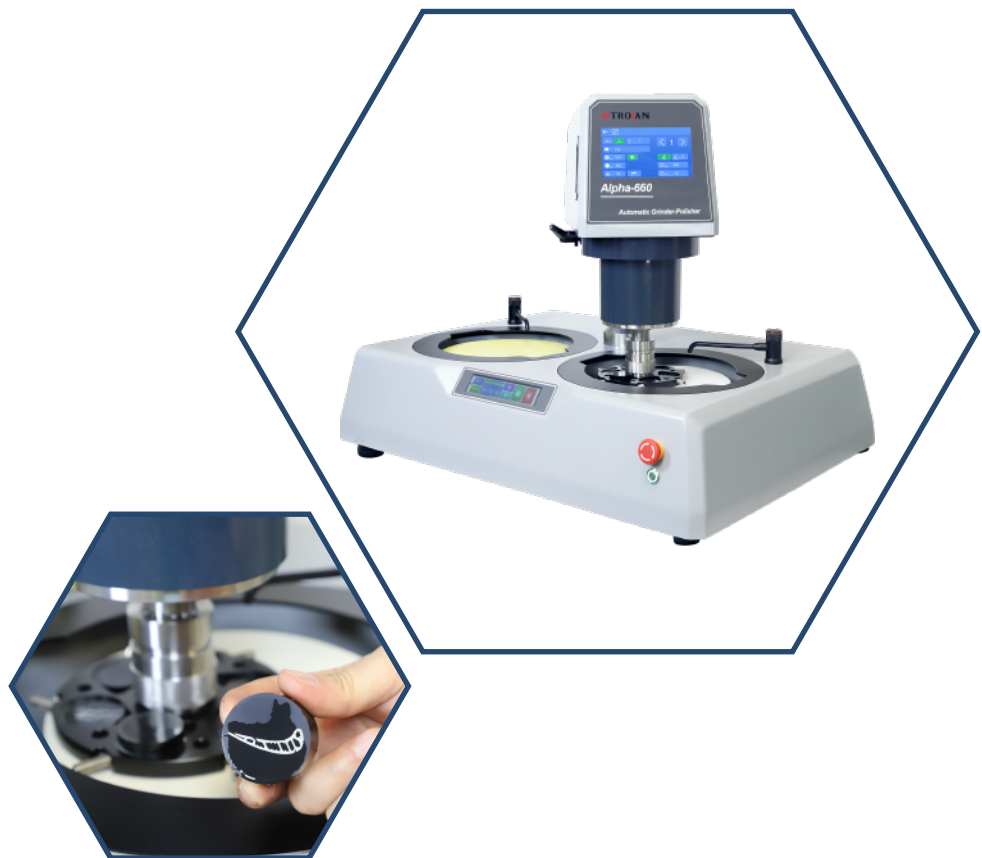

2025
03 / 04

TROJAN (SUZHOU) MATERIAL TECHNOLOGY CO.,LTD

Alpha-660 Dual-Disc Automatic Grinder and Polisher

The Alpha-660 is a dual-disc manual-automatic grinder and polisher designed for precision sample preparation. It features independent motor drives for the head and polishing discs, supporting disc diameters from 200mm to 300mm. The machine offers two pressure modes: central and individual loading. An optional peristaltic pump system with up to four channels is available for precise dispensing of polishing solutions. The system allows for the programming and storage of up to 20 commonly used grinding and polishing procedures, ensuring consistent and reproducible results with ease and efficiency.





Features:

- ◆ Constructed with a cast aluminum shell and SMC standard loading cylinders, the design is more reliable, robust, and durable;
- ◆ Features a 7-inch touchscreen panel with a user-friendly and intuitive interface;
- ◆ Dual-touchscreen design with independent control of dual motors, allowing simultaneous automatic and manual grinding;
- ◆ The sample holder height is adjusted using a dual-hand safety button, requiring simultaneous activation;
- ◆ Offers both central and individual loading options; central loading allows for precise, quantitative grinding;
- ◆ Optional peristaltic pump system with four channels for precise solution dispensing;
- ◆ Includes 20 pre-programmed grinding and polishing processes, with support for custom parameter storage;
- ◆ Equipped with a large-diameter side-outlet drainage system to prevent clogging and facilitate easy drainage;
- ◆ Features a detachable faucet design for easy cleaning and maintenance of the discs
- ◆ Built-in automatic temperature-sensitive cooling fan reduces failure rates and extends service life.



Dimension Figure


Unit: mm

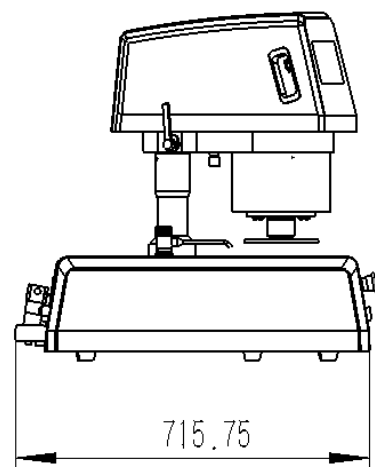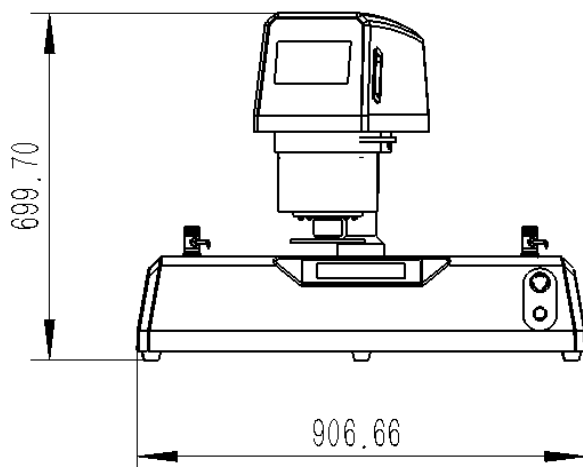

Technical Data:


Model		Alpha-660
Grinding Head	Loading force	Individual loading & Central loading
	Central loading	20-360N
	Individual loading	5-60N
	Power	220V, 750W
	Speed	20-150rpm
	Direction	CCW/CW
Main Unit	Power	220V, 2x750W
	Speed	50-600rpm
	Diameter	10"/254 (8"/203mm, 12"/305mm optional)
	Direction	CCW
Quantitative Function	Accuracy	±0.1mm
	Quantitative Height	0.1mm-10mm
Dosing System (optional)	Type	Peristaltic pump, DC24V, 10W
	Flow	Maximum flow rate: 52ml/minute
	Amount	4
Clamping Disc	Individual loading	Standard: 6xΦ30mm(6xΦ32mm,6xΦ40mm,6xΦ50mm optional)
	Center loading	Standard: 6xΦ32mm(6xΦ25mm,6xΦ40mm,6xΦ50mm optional)
Environment	Temperature	-10°C - 60°C
	Humidity	0-95%RH
Dimension	WxDxH	716x907x700mm
Weight		100kg

Control Panel Introduction:

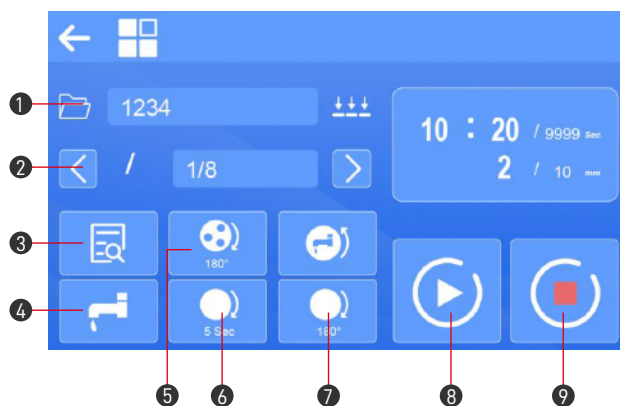

Automatic mode:


- ① Select program
- ② Select subroutine
- ③ Preview program
- ④ Turn on or off the faucet manually
- ⑤ Place the sample by manually rotating the clamping disk by 180°
- ⑥ Manual dumping
- ⑦ Rotate the working platen by 180° for sample removal
- ⑧ Start
- ⑨ Stop

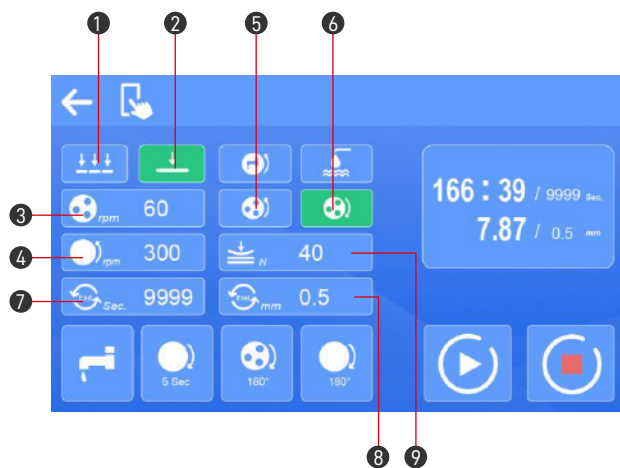

Manual mode:

- ① Individual loading
- ② Central loading
- ③ Speed setting for working platen (50-600) RPM
- ④ Speed setting for clamping disk (20-150) RPM
- ⑤ Forward rotation
- ⑥ Reverse rotation
- ⑦ Setting for running time (0-9999) sec
- ⑧ Quantitative grinding value (0.1-10) mm (optional)
- ⑨ Pressure setting, central pressure: (20-360) N. Individual pressure: (5-60) N

Related Consumables & Accessories:

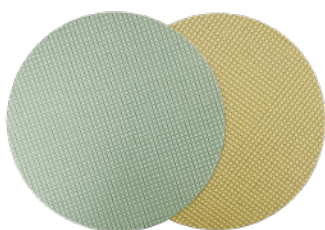

DiaRe Diamond Disc



Abrasive Paper



Polishing Cloth



Diamond Suspension



Alumina/Silica suspension



Adapter Plate System

Standard Accessories:

1	Working platen		10" (254mm)	2pcs
2	Individual clamping disc and connector		Φ30 Specimen holding platen	1set
3	Inlet pipe		Φ8 PU pipe + filter	2pcs
4	Outlet pipe		Corrugated Tube with an Inner Diameter of Φ32mm and Clamp	2pcs
5	Intake Pipe		Φ6mm PU pipe	1pc
6	Abrasive paper		180#(250mm),400#(250mm), 800#(250mm),1200#(250mm), 2500#(250mm)	10pcs each
7	Magnetic disk		Magnetic disk (250mm)	2pcs
8	Conversion disk		AD-PB (250mm)	2pcs
9	Polishing cloth		ET-JP (250mm) YS-JP (250mm)	2pcs each
10	Power cord		1.5M	1pc

Optional Accessories:

1	Working platen		10" (254mm) or 12" (305mm)	
2	Magnetic disk		Magnetic disk (250mm or 300mm)	
3	Conversion disk		AD-PB (250mm or 300mm)	
4	Individual clamping disc and connector		6xΦ30mm, 6xΦ32mm, 6xΦ40mm, 6xΦ50mm, or other customized	
5	Teardrop-shaped clamping disc with leveling base		6xΦ30mm, 6xΦ32mm, 8xΦ40mm, 6xΦ50mm, or other customized	

Optional Accessories:

6	Specimen Leveling Platen		Φ167mm
7	FLAT Flattening		FLAT-1
8	ADS-4 Dispenser		500mlx4, flow rate :8-52ml/minute
9	Air Compressor		220V, 1.1kW, Exhaust Pressure 0.8MPa, Capacity 40L

Purchase list:

DiaNi Diamond Disc:

	8in DiaNi (200mm)	10in DiaNi (250mm)	12in DiaNi (300mm)	Package
80	04.21.110S	04.21.210S	04.21.310S	1PC
120	04.21.120S	04.21.220S	04.21.320S	1PC
200	04.21.125S	04.21.225S	04.21.325S	1PC
400	04.21.130S	04.21.230S	04.21.330S	1PC
800	04.21.140S	04.21.240S	04.21.340S	1PC
1200	04.21.150S	04.21.250S	04.21.350S	1PC

DiaRe Diamond Disc:

	8in DiaRe (200mm)	10in DiaRe (250mm)	12in DiaRe (300mm)	Package
80	04.22.110S	04.22.210S	04.22.310S	1PC
120	04.22.120S	04.22.220S	04.22.320S	1PC
200	04.22.125S	04.22.225S	04.22.325S	1PC
400	04.22.130S	04.22.230S	04.22.330S	1PC
800	04.22.140S	04.22.240S	04.22.340S	1PC
1500	04.22.150S	04.22.250S	04.22.350S	1PC
3000	04.22.160S	04.22.260S	04.22.360S	1PC
5000	04.22.170S	04.22.270S	04.22.370S	1PC

Silicon Carbide Disc:

	8in PSA (200mm)	10in PSA (250mm)	12in PSA (300mm)	8in Plain Back (200mm)	10in Plain Back (250mm)	12in Plain Back (300mm)	Package
120	03.04.203	03.04.503	03.04.803	03.03.103	03.03.403	03.03.703	50PCS
180	03.04.204	03.04.504	03.04.804	03.03.104	03.03.404	03.03.704	100PCS
400	03.04.207	03.04.507	03.04.807	03.03.107	03.03.407	03.03.707	100PCS
800	03.04.209	03.04.509	03.04.809	03.03.109	03.03.409	03.03.709	100PCS
1200	03.04.211	03.04.511	03.04.811	03.03.111	03.03.411	03.03.711	100PCS
2000	03.04.213	03.04.513	03.04.813	03.03.113	03.03.413	03.03.713	100PCS
2500	03.04.214	03.04.514	03.04.814	03.03.114	03.03.414	03.03.714	100PCS
2400	03.04.215	03.04.515	03.04.815	03.03.115	03.03.415	03.03.715	100PCS
4000	03.04.216	03.04.516	03.04.816	03.03.116	03.03.416	03.03.716	100PCS

Purchase list:

Alumina Disc:

	8in PSA (200mm)	10in PSA (250mm)	12in PSA (300mm)	8in Plain Back (200mm)	10in Plain Back (250mm)	12in Plain Back (300mm)	Package
60	03.02.401	03.02.501	03.02.801	03.01.101	03.01.101	03.01.701	50PCS
120	03.02.403	03.02.503	03.02.803	03.01.103	03.01.103	03.01.703	100PCS
180	03.02.404	03.02.504	03.02.804	03.01.104	03.01.104	03.01.704	100PCS
240	03.02.404-1	03.02.504-1	03.02.804-1	03.01.104-1	03.01.104-1	03.01.704-1	100PCS
320	03.02.405	03.02.505	03.02.805	03.01.105	03.01.105	03.01.705	100PCS
400	03.02.406	03.02.506	03.02.806	03.01.106	03.01.106	03.01.706	100PCS
800	03.02.408	03.02.508	03.02.808	03.01.108	03.01.108	03.01.708	100PCS
1200	03.02.410	03.02.510	03.02.810	03.01.110	03.01.110	03.01.710	100PCS

Adapter Plate System:

	AD-PSA Plate	AD-Foil Plate	AD-PB Plate	Magnetic Plate	Package
200mm	04.08.410	04.08.422	04.08.310	04.08.010	1PC
250mm	04.08.413	04.08.430	04.08.320	04.08.020	1PC
300mm	04.08.420	04.08.432	04.08.330	04.08.030	1PC

Polishing Cloth:

		8in (203mm)	10in (254mm)	12in (305mm)	Package
Rough Polishing	YS-JP White Fiber Polishing Cloth(PSA)	05.07.710	05.07.720	05.07.730	10pcs
	GF-JP White Synthetin Silk Polishing Cloth (PSA)	05.10.010	05.10.020	05.10.030	10pcs
Fine Polishing	SR-ZP Red Synthetic Fabric Polishing Cloth (PSA)	05.02.150	05.02.160	05.02.170	10pcs
	SC-JP White Natural Silk Polishing Cloth (PSA)	05.05.510	05.05.520	05.05.530	10pcs
	CS-JP Blue Synthetin Silk Polishing Cloth (PSA)	05.12.110	05.12.120	05.12.130	10pcs
Final Polishing	ZN-ZP Black Polyurethane Polishing Cloth (PSA)	05.01.110	05.01.120	05.01.130	10pcs
	ET-JP White Synthetic Velvet Polishing Cloth (PSA)	05.03.210	05.03.220	05.03.230	10pcs
	YR-JP Brown Synthetic Rayon Polishing Cloth (PSA)	05.04.310	05.04.320	05.04.330	10pcs
	Final Red Final Precision Polishing Cloth(Adhesive Back)	05.13.110	05.13.120	05.13.130	10pcs

Diamond polishing suspension and paste:

Diamond Suspension		Micrometer	Diamond Paste		Diamond Spray
Monocrystalline (500ml)	Polycrystalline (500ml)		Monocrystalline (5g)	Polycrystalline (10g)	Monocrystalline (350ml)
06.01.110-1	06.02.210	0.25um	—	—	—
06.01.120-1	06.02.220	0.5um	—	—	06.03.115
06.01.130-1	06.02.230	1um	06.10.030(Blue)	06.10.110(Black)	06.03.320
—	—	2.5um	—	—	06.03.325
06.01.140-1	06.02.240	3um	06.10.040(Green)	06.10.130(Black)	—
—	—	3.5um	—	—	06.03.330
—	—	5um	—	—	06.03.335
06.01.150-1	06.02.250	6um	06.10.050(Yellow)	06.10.160(Black)	—
06.01.160-1	06.02.260	9um	06.10.060(Red)	06.10.190(Black)	—
—	—	10um	—	—	06.03.340

Purchase list:

Precision polishing suspension and polishing powder:

	0.05μm	0.3μm	1μm
Alumina Suspension(500ml)	06.04.210	06.04.220	06.04.230
Alumina Polishing Powder(1LB)	06.13.011	06.13.021	06.13.031

	500ml	1000ml	5000ml
SO T401 Non-Crystallizing Silica(50nm)	06.05.310	06.05.320	06.05.330
100C30 silica suspension(50nm)	06.08.310	06.08.311	06.08.312
SO A539 silica suspension(50nm)	06.07.520	06.07.523	

Polishing lubricants:

	5000ml	500ml	1000ml
PL-W Polishing Lubricants	-	06.09.510	06.09.520
PLW-B Polishing Lubricants	06.09.630	06.09.640	06.09.650

Grinding and Polishing Disc:

	200mm	250mm	300mm	350mm
POS Grinding and Polishing Disc(1pcs)	04.09.110	04.09.115	04.09.120	04.09.125
POH Grinding and Polishing Disc(1pcs)	04.09.130	04.09.135	04.09.140	04.09.145



TROJAN (SUZHOU) MATERIAL TECHNOLOGY CO.,LTD

Tel: +86-512-65283666

Email: sales@trojanchina.com

Web: www.trojanworld.com

Add: Building 26, No.228 Gongtang Road, Luzhi Town, Wuzhong Suzhou, China

